
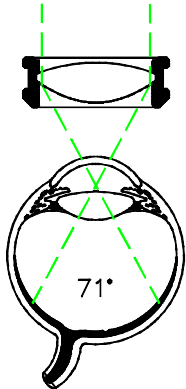


Ocular Indirect Vitrectomy Lenses

 CE	Product Code	Image Mag	Static FOV	Dynamic FOV	OIV-60	OIV-132
		OIV-132	.45x	99°	135°	

Design

Designed to be used on the Ocular Landers Wide Angle Surgical Viewing System (OSVS) in conjunction with an Inverter Vitrectomy System.

Attaches to the OSVS with the Ocular 132D Indirect Vitrectomy Lens Holder (OIV-H132)

Non-contact design allows patient's eye to be rotated freely to view the peripheral retina and vitreous

Steam sterilizable; sterilizing case included

Designed to allow a clear view in the fluid or air filled eye

Cleaning and Disinfection

See Cleaning Method 3

